

Stools and Chairs



Ritter[®]
by MIDMARK

Let our seating
products provide
the comfort and
support you
deserve.

You are up and down all day, every day - so what you need most in seating is comfort and support. Our seating products offer wide, amply padded seats for greater comfort while a strong base structure provides superb stability.

The Ritter brand has been providing quality products to the healthcare industry for over 50 years and is recognized as an industry leader...our seating products are durable and built to last.

Coordinating your rooms with Ritter and Midmark tables and casework, will provide a warm, inviting look that your staff and patients will recognize and appreciate.



Stools & Chairs



Designed and manufactured to exceed your expectations...

It begins with a wide, five-leg caster base that provides support and stability. The star-shaped base allows closer accessibility to the patient while heavy duty, dual-wheel, hooded casters enable stools to maneuver quietly over all types of surfaces.

Adding to the strength and stability of our stools is the solid steel machine screw on manual models, or air cylinder on pneumatic models. The air cylinder further acts as a built-in shock absorber, taking the abuse of a hard day's work instead of you, and presents a simple and efficient method for adjusting the seat height; ideal for multi-physician practices.

Comfort is provided in the form of a wide, amply padded seat section constructed of dense foam. A beadless upholstery seam eliminates the seat edge from digging into your legs. Add a backrest to promote better posture and support to leave you feeling less fatigued at day's end.

With many choices, colors and accessories, there is certain to be a stool just right for you!

276 Air Lift Stool

Our best selling stool! The 276 has become the industry standard in healthcare facilities.



Classic Series Stools

The healthcare industry's most demanded stools, and the foundation of our seating product line. The Classic Series features a polished aluminum/black accented five-leg caster base and hood, and a hand ring or foot lever height control on pneumatic models. The seat height of the 274 (without back) and 275 (with back) adjusts manually while the 276 (without back) and 277 (with back) adjusts pneumatically. The 278 (foot lever) and 279 (hand ring) offer pneumatic height adjustment and a surgical armrest to stabilize the surgeon's arm during delicate procedures.



Stools

427 Air Lift Physician Stool

Same features as the 425, but with a foot lever instead of a hand ring. A touch of your heel to the foot lever pneumatically adjusts the height, allowing you to maintain a sterile field. The 427 is shown with optional back kit that adjusts in, out, up and down; and optional arm kit that adjusts in, out and angularly.

** 427 shown with optional back and arm kit accessories*



425 Air Lift Physician Stool

The ergonomic design of the 425 features a contoured seat that molds to your body shape, making it the most comfortable stool you will ever sit on! A simple lift of the hand ring pneumatically adjusts the seat height, allowing you to find your most comfortable working height effortlessly and efficiently. Plastic covers hide the steel substructures and five-leg caster base to provide a warm aesthetic appearance. The 425 and 427 redefine the expectations of healthcare professionals.



Value Series Stools

Economy version of the popular Classic Series Stools without compromising reliability and style. The Value Series offers a black, plastic composite five-leg caster base and hood, and a single lever height control on pneumatic models. The seat height of the 270 (without back) and 271 (with back) adjusts manually while the 272 (without back) and 273 (with back) adjusts pneumatically.



680 Side Chair

Soft and comfortable seating for your patients in the exam room, waiting room and front office. The 680 comes with all the extras: contoured seat and back section for added support; sturdy, powder coated black steel tube frame with leveling feet for increased stability and a “wall-saver” kickback leg design that allows the chair to be close to the wall without damaging it. Available with or without high-impact plastic chair arms.

281 Blood Drawing Chair

The slightly angled back and adjustable chair arms offer patient comfort and security. Arms are equipped with quick release mechanisms for simple and efficient positioning of their height and inward/outward angle. The 281 is available with a drawer, ideal for storing supplies and providing additional work surface area.



Stools and Chairs

S P E C I F I C A T I O N S

Classic Series Stools

Seat Height:
18" - 24.5" (274-278)
19" - 25.5" (279)

Seat Diameter: 16.25"

Seat Cushion Thickness: 3.25"

Back Cushion Dimensions:
14" W x 7.5" H (275, 277)

Caster Base Diameter: 21.5"

Weight Capacity:
350 lbs. (274, 276)
225 lbs. (275, 277, 278, 279)

Value Series Stools

Seat Height:
18" - 24.5" (270-273)

Seat Diameter: 16.25"

Seat Cushion Thickness: 3.25"

Back Cushion Dimensions:
14" W x 7.5" H (271, 273)

Caster Base Diameter: 23"

Weight Capacity:
300 lbs. (270, 272)
225 lbs. (271, 273)

Classic and Value Series Options

Locking casters
Soft rubber casters
Glides
6" Lab height extension
Foot ring

425/427 Air Lift Physician Stools

Seat Height:
17.5" - 24" (425)
17" - 23" (427)

Seat Cushion Dimensions:
16.5" W x 17" D

Seat Cushion Thickness: 3"

Back Cushion Dimensions:
13" W x 9.5" H
(with back kit accessory)

Caster Base Diameter: 23"

Weight Capacity: 225 lbs.

281 Blood Drawing Chair

Seat Height: 20"

Seat Cushion Dimensions:
20" W x 16.5" D

Seat Cushion Thickness: 2"

Back Cushion Dimensions:
20" W x 16.5" H

Inside Arm Width: 21"

Outside Arm Width: 32"

Storage Drawer Capacity:
8.5" W x 15.25" D x 4.5" H

Weight Capacity: 400 lbs.

680 Side Chair

Seat Height: 17"

Seat Cushion Dimensions:
18" W x 19.5" D

Seat Cushion Thickness: 2"

Back Cushion Dimensions:
18" W x 13.5" H

Inside Arm Width: 17.5"
(on model with arms)

Outside Arm Width: 23.5"
(on model with arms)

Weight Capacity: 400 lbs.


Upholstery Colors:



Ritter: Navy, Dusty Blue, Moss, Clay, Pebble Grey, Shadow, Perfect Plum and Black



Midmark Premium Solid:
Blueberry, Silver Sage, Tea Green, Fossil Grey, Cashmere Blue, Terra Cotta and Black

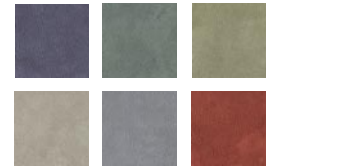

Midmark Premium Pattern:
Belagio, Siena, Tuscany, Firenze, Capri and Treviso



Midmark Ultra-Premium:
Ultra-Pacific Blue, Ultra-Spruce, Ultra-Olivine, Ultra-Arctic, Ultra-Sky, Ultra-Canyon and Ultra-Raven Wing

*There will be an additional charge for Midmark Ultra-Premium.

Midmark is an ISO 9001 Certified Company.

For more information or a demonstration, contact your Ritter or Midmark Dealer, or contact Midmark Corporation, Versailles, Ohio.
Call: 1-800-MIDMARK
Fax: 1-800-365-8631
midmark.com



Classic and Value Series Backrest



Classic Series Glides



Value Series Glides



Classic and Value Series Locking Casters



Classic Series Soft Rubber Casters



Value Series Soft Rubber Casters



Classic Series 277 shown with optional 6" lab height extension and foot ring.

Ritter[®]
by MIDMARK

© Midmark Corporation 2005

Products subject to improvement changes without notice.
Litho in U.S.A. 007-0558-00 Rev. C1 (3/07)