


Perfect for large surgery suites that use more than one type of anesthetic agent

- Three gas machine (O₂/N₂O/Air)
- Holds two vaporizers
- Large work area with a writing table and footrest
- Includes a recessing absorber

► A Quick Look

- Mechanical Anti-Hypoxic Device and Air/N₂O interlock
- Ventilator options (not included): DRE AV800, Ohmeda 7000 or 7800
- Selectatec®-compatible vaporizer mounting system
- Holds up to two vaporizers
- CO₂ absorber slides under work shelf
- Double-canister CO₂ absorber standard
- Specifically designed for low flow
- Left or right handed layout
- Contemporary ergonomic design
- Recessing absorber and ventilator
- Only two preventive maintenances required each year

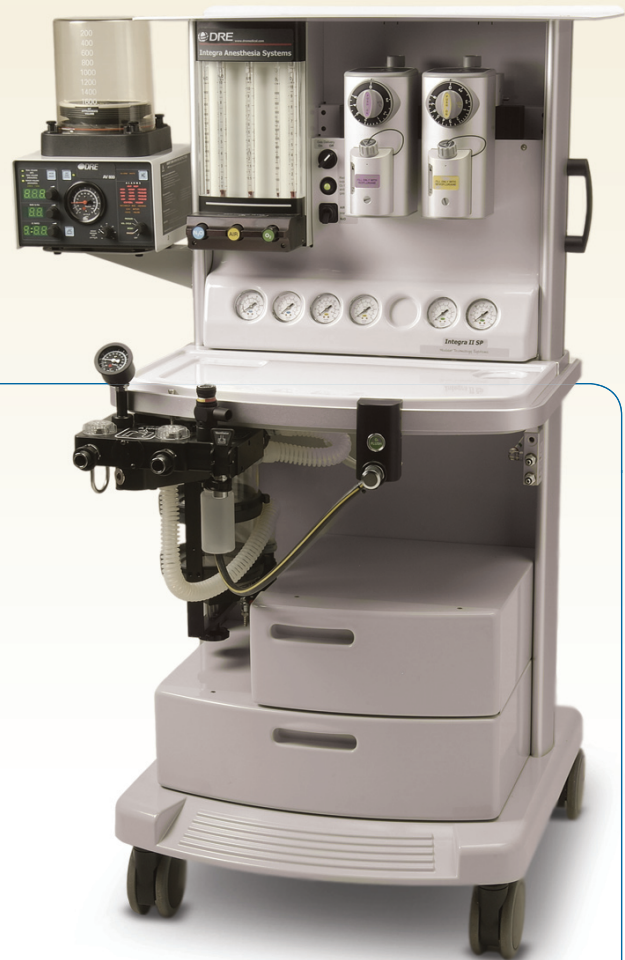

Physical Dimensions

Height: 54" (137 cm)
Width: 26" (66 cm)
Depth: 26" (66 cm)
Weight (with all options): 200 lbs. (75 kg)

Power

Power Requirements: 110/220 VAC
Auxiliary Electric Outlets: (4) 13 Amp

Display

Top Shelf: 24" x 16" (61 cm x 40 cm)
Work Surface: 24" x 14.5" (61 cm x 37 cm)
Drawers: 6" x 15" x 17"/max: 3
(15 cm x 38 cm x 43 cm)
Writing Tablet: A4 size

Order now

to receive a **two-year manufacturer's warranty**
and a **two-year parts warranty**