



High quality function elements designed with practical, dependable materials and craftsmanship.


The DRE Lucerne is the effective, economical choice for your surgery suite

Features

- The electro-hydraulic compact mechanism guarantees high operation safety and precise patient positioning.
- The radiolucent table top permits maximum full body use for a C-arm X-ray unit.
- The table top is divided into five sections with a complete detachable mattress to provide high performance flexibility. The head and leg sections can also be removed.
- The remote control—with coiled cable—adjusts the table height, lateral tilt, back and pelvis sections.
- The maneuverable and stable base has moveable castors which can be rolled both longitudinally and laterally.
- Equipped with the optional ups battery box that supplies energy for up to 100 operations.
- An electro-hydraulic system is integrated into the table base. This control system is specially designed with an overload protection device that can prevent equipment damage by abnormal motor pump overload.

► Specifications

ElecAC input: 110V or 220V 50/60Hz
Remote controller power: DC24V
Battery: 50AH, N50 (2PS)
Table weight: 375 lbs
Mattress thick: 2.2 in
Patient weight capacity: 397 lbs
Table height adjustment (w/o mattress): min. 28 in, max. 40 in
Trendelenburg: 35°/35°
Fowler: 15°/60°
Lateral tilt adjustment:
Left: 20°, **Right:** 20°
Back section adjustment:
Up: 60°, **Down:** 22°
Head Plate Adjustment:
Up: 20°, **Down:** 90°
Pelvis section adjustment:
Up: 40°, **Down:** 20°
Leg plate adjustment:
Up: 20°, **Down:** 90°

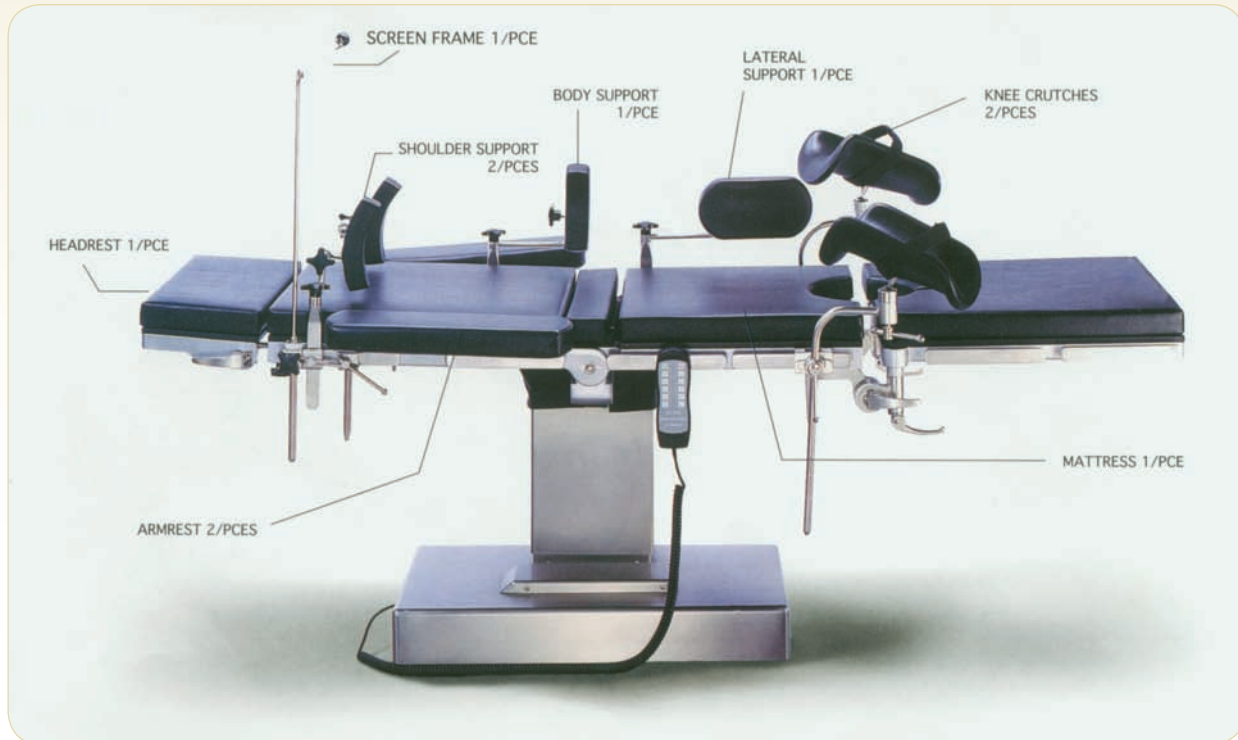

Remote Control

- Automatic resetting of the table
- Back section tilt
- Pelvis section tilt
- L.E.D. indicating ON/OFF
- Table raising or lowering
- Lateral tilt
- Electro hydraulic locking device



Radiography Table Top

Use for C Arm X-Ray



Additional features

- Emergency safety button
- UPS battery plug
- Fuse
- Power input
- ON/OFF switch