


LIFT CAPACITY: 1,000 POUNDS

Strength and dexterity with a full-length radiolucent top.
Sleek, easy-to-clean stainless steel construction.


- » Powered surgical table with a 180° rotating top — provides reliable C-arm access.
- » Full-length radiolucent top and superior imaging from an offset cylindrical column.
- » Ample lift capacity: 1,000 pounds.
- » Smooth steering castors provide enhanced mobility.
- » Additional features and options: Battery back-up, removable head and leg section, 2" pressure-relieving mattress, and auto-leveling brake system.

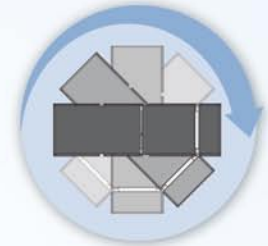

180° Rotating Top
for precise patient positioning




12" horizontal slide

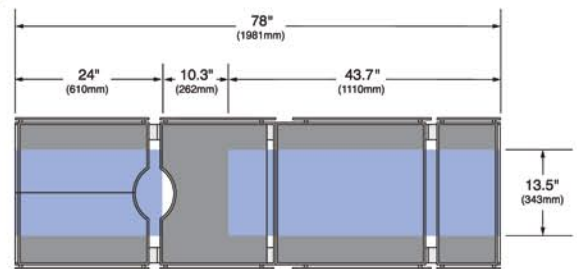

X-ray via cassette tunnel, and C-arm imaging

Table Traversed Toward Head →

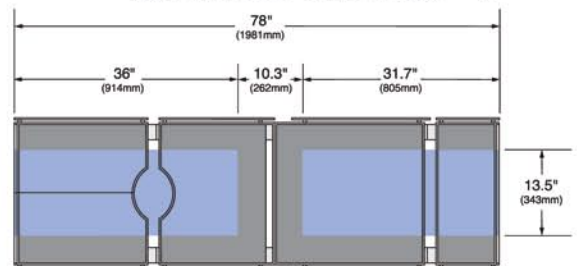

X-ray via cassette tunnel, and C-arm imaging

← **Table Traversed Toward Foot**

The 12" horizontal sliding top offers superior full-body imaging access and complete surgical control.



Easy to position using hand control or optional foot control

DRE Paramount 1000SL

Surgical Table

Equipment for the way *you* operate

Specifications

Tabletop

Length	78"
Width	19.75"
Max. height	43"
Min. height	29"

Weight Information

Lift Capacity	1,000 lbs.
Articulation	800 lbs.
Flex	600 lbs.




Maximum allowable weight of:

Back section	150 lbs.
Leg section	100 lbs.
Head section	28 lbs.

Positioning

Lateral tilt	30°
Trendelenburg	30°
Back section	+90°, -30°
Head section	+30°, -90°
Leg section	+0°, -105°
Table rotation	180°
Self-leveling brakes	Powered

Specifications subject to change without notice.

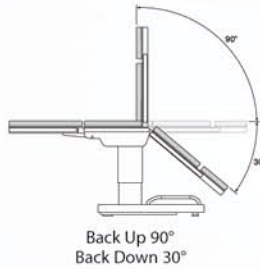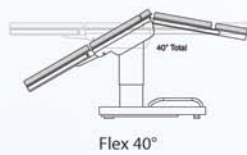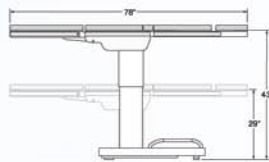

Neurological Seated Position



Nephrectomy Procedure



Laminectomy Procedure



Gynecology and Urology Procedures

Traction Assembly Option

Requires only one person for installation.



For orthopedic procedures:

- Fibia and Tibia • Lateral • Hip Fractures

Shoulder Chair Option

Attaches effortlessly to the table



Lateral support and removable wings for access

Accessories subject to availability.