


High Beam TriLens



NEVER COMPROMISE
START WITH YOUR VISION

We raised the bar for surgical headlights

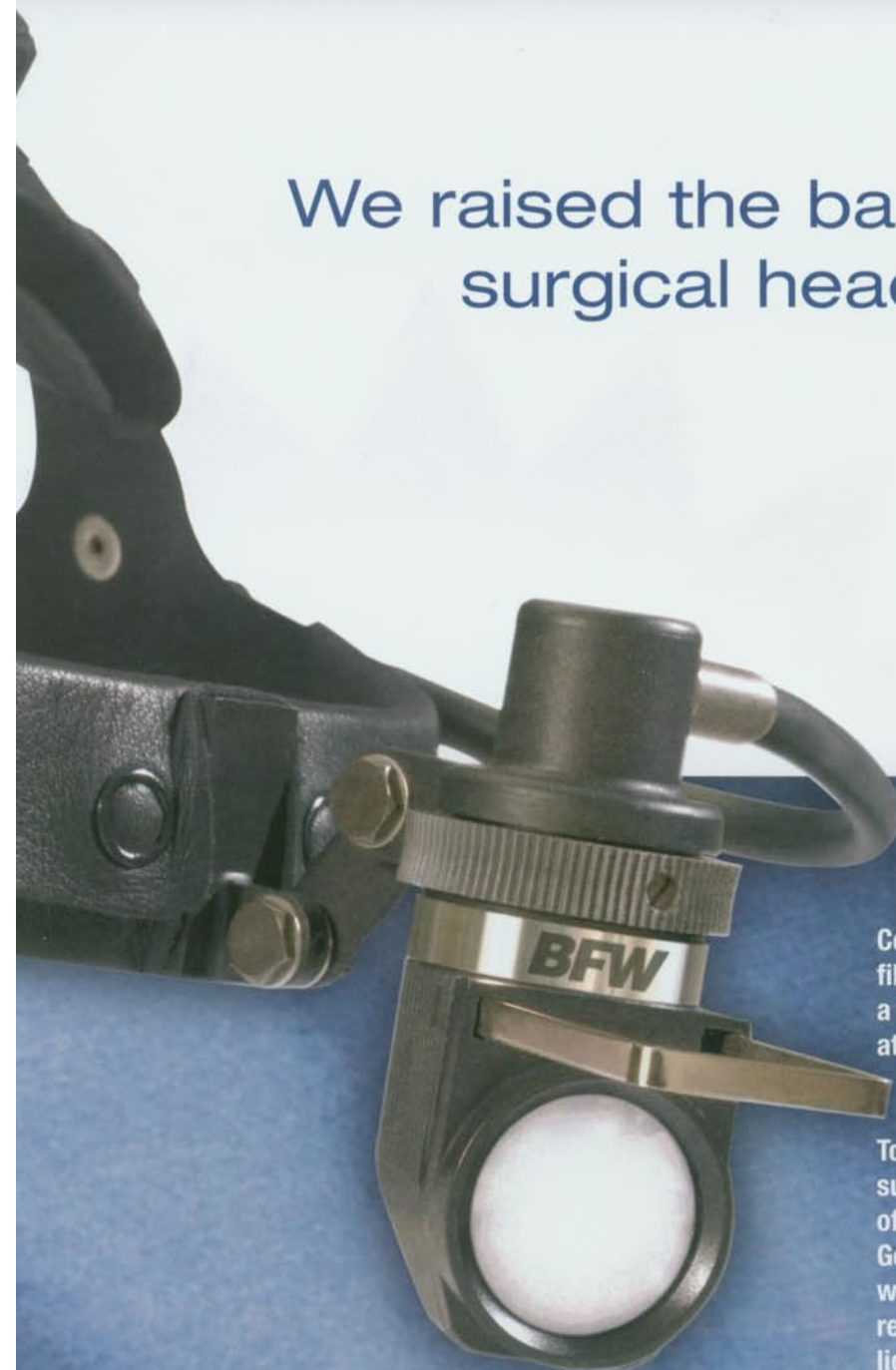

Comfort is a given with our High Beam TriLens fiber optic headlight. Incomparably, it provides a clear and remarkably brilliant 11.5 cm spot at a 40cm working distance – the largest in its class.

To achieve this unequaled performance in a surgical headlight, we craft our optics with three of the finest Mag-Fluoride coated lenses and a German-engineered iris. We mold the optic crown with DuPont Minlon®, one of the preeminent heat resistant materials known today. Minlon® is lightweight and ensures safety and durability when using xenon light.

Our 4.75 mm fiber optic cable holds 9,000 glass strands sheathed in our new ultra light BFW™ BLUE Opaque sheath, affording a lower profile than cables with half the light delivering capacity.

Choose our form-fitted, sturdy and dual-adjustable Classic headband or our soft Breath-O-Prene® Sport headband. Both offer unmatched comfort, performance, safety and durability.

Unequaled performance
in a surgical headlight



1-800-477-2006 • www.dremed.com

High Beam TriLens



NEVER COMPROMISE
START WITH YOUR VISION



High Beam Classic

The High Beam TriLens Fiber Optic headlight performs flawlessly with our Maxenon™ Xi300 and Maxenon™ PowerPlus xenon light sources and with our Maxillume™ or MiniMax™ T3 Quartz Halogen series.

Specifications

Product Codes

- MID 7800
High Beam TriLens Classic with 7' Cable
- MID 7810
High Beam TriLens Classic with 10' Cable
- MID 7830
High Beam TriLens Sport with 7' Cable
- MID 7835
High Beam TriLens Sport with 10' Cable

Accessories

- MID 7704
Positioning Bar
- MID 8600
7' Replacement Cable
- MID 8610
10' Replacement Cable
- MID 7825
Monocular Optic (Lawson Lite)
- MID 4450
Headlight Carrying Case



High Beam Sport



1-800-477-2006 • www.dremed.com